


ragmask v2.2b1

ragmask.com


1. Print just this page.

You don't need the other pages, save your paper.

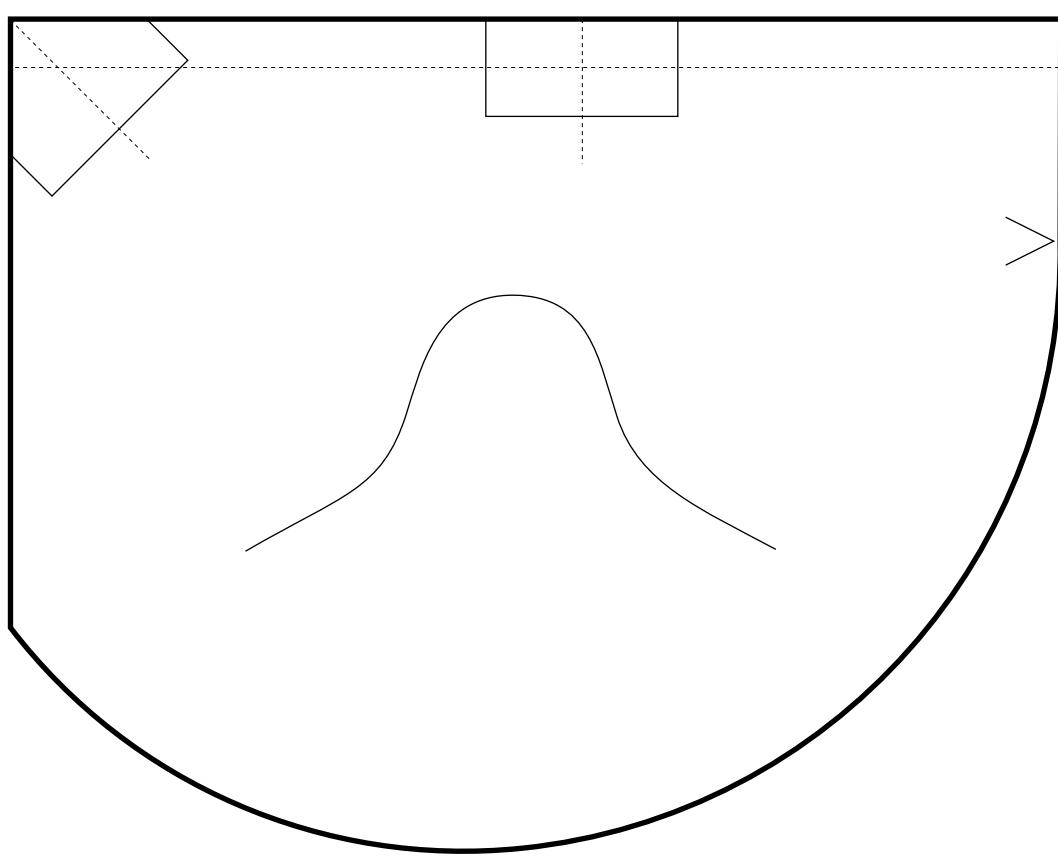


2. Check the scale.

Make sure it printed at the right size. A real ruler should match the ruler below.

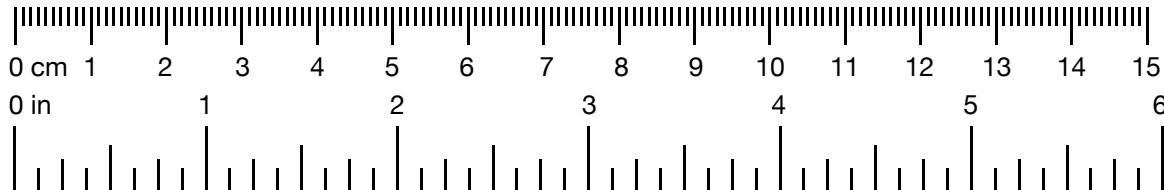

3. Cut materials.

11"x11" (28x28cm) fabric, 1 piece per mask
1"x16" (2.5x40cm) fabric for straps, 4 pieces per mask

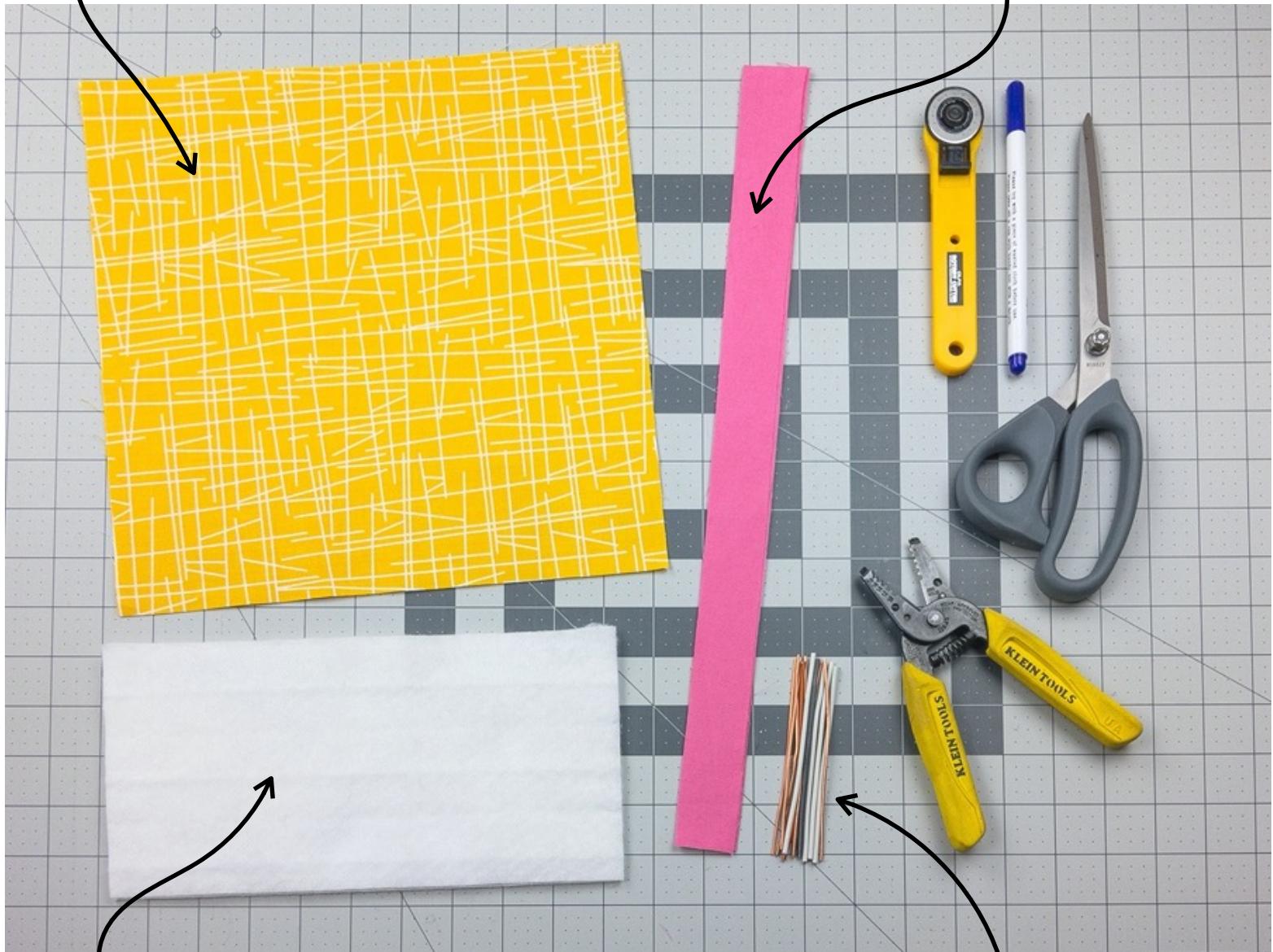


10"x5" (26x13cm) filter (optional)*, 1 piece per mask
4" (10cm) metal wire, 1 piece per mask, 14ga copper wire works well

* Effectiveness of various filter materials forthcoming.

See ragmask.com/materials for latest info.

Even no filter (just fabric) can be more effective than no mask at all.

4. Fold fabric over filter.

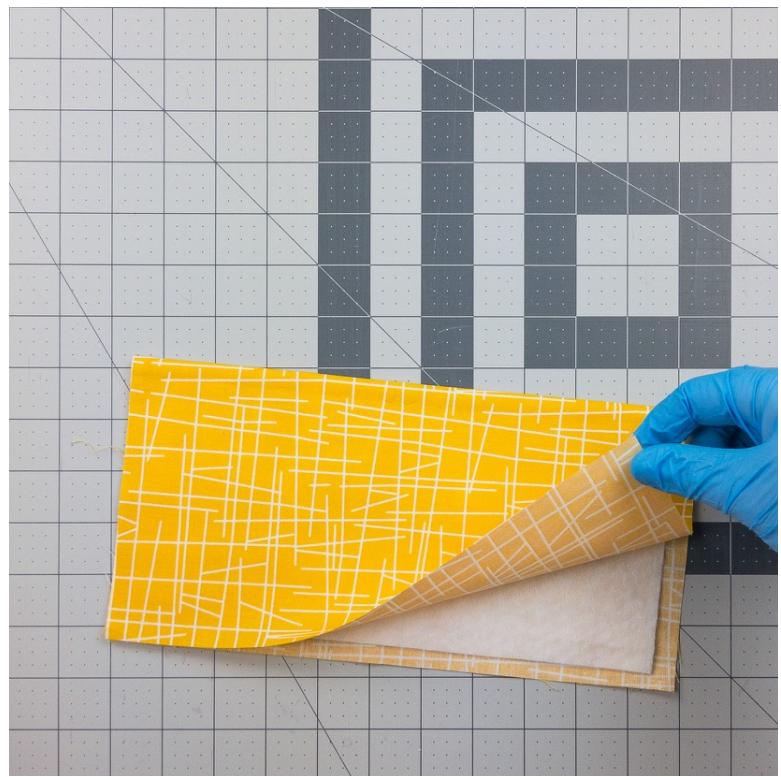


5. Sew around 3 sides.


~3/16" (0.5cm) seam allowance, quick > pretty.



6. Fold.

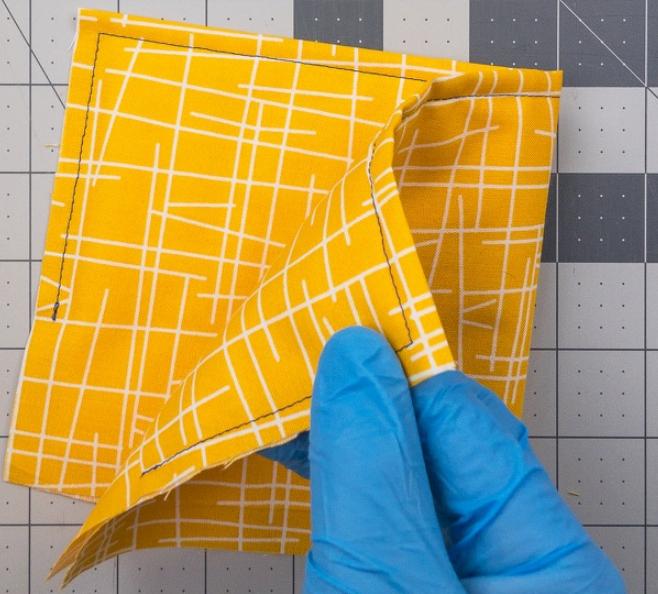

7. Mark.

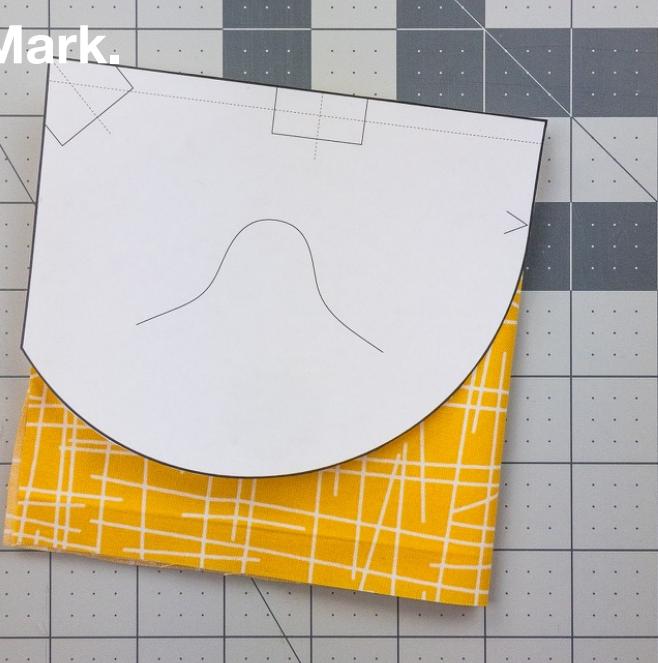

8. Sew.




9. Cut.



10. Flip inside out.



11. Mark strap locations.

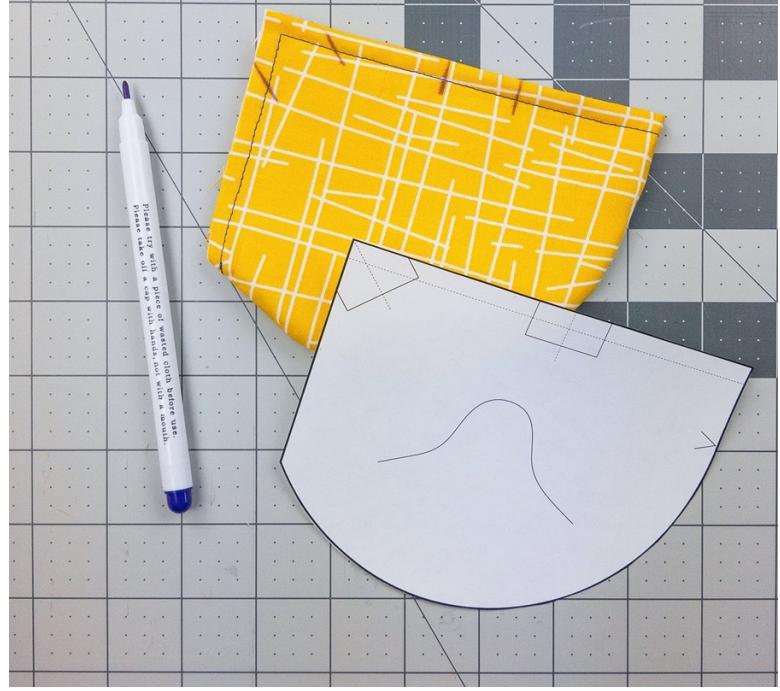


12. Sew straps on outside face of mask. See last page for an even faster way to install straps.





13. Insert wire and center over nose.

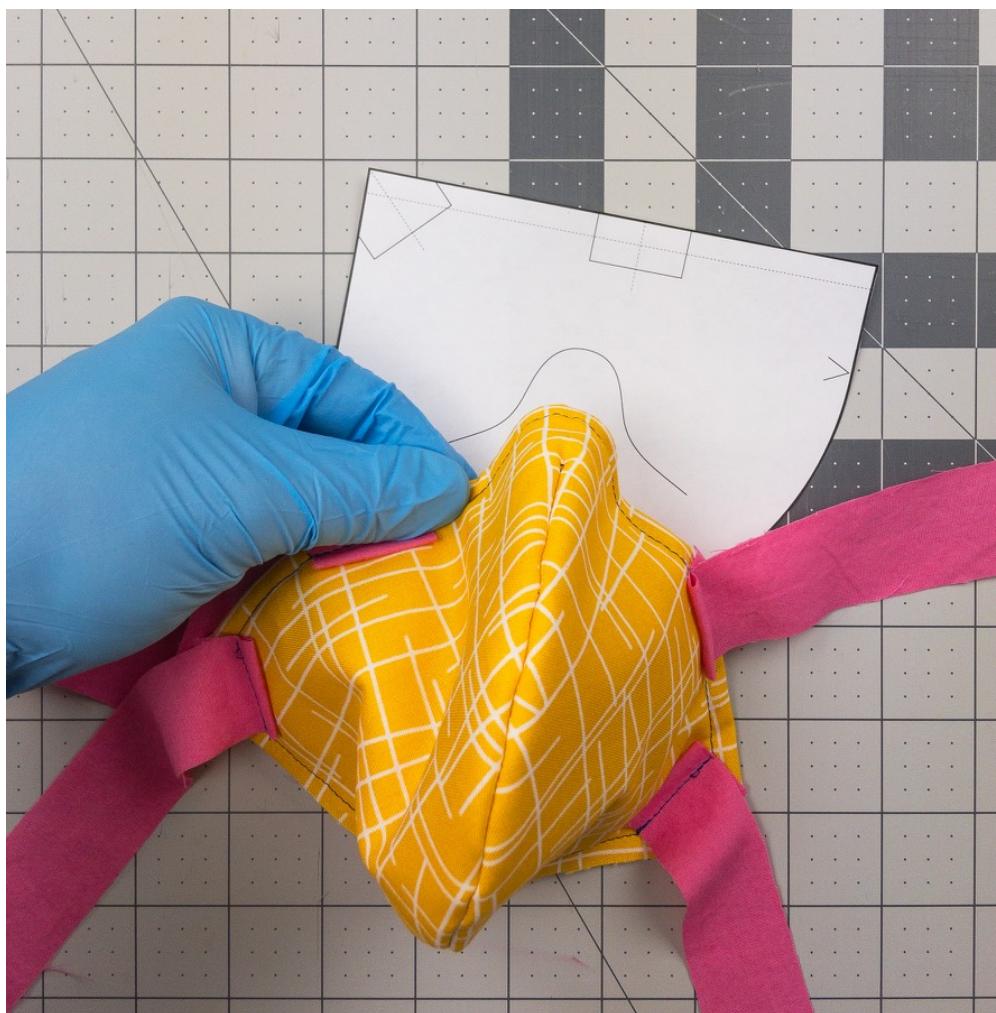

14. Form wire into approximate nose shape.

Wear

ragmask.com



Use and care details forthcoming...

Pro Tip: Fast Straps

More tips at ragmask.com/pro



1. Print another copy of Page 1, trim template. After making first fold, mark centers of straps.



2. Sew straps on as part of the first stitch. (Backstitch over straps for added durability).



3. Fold in half (straps on the inside of the fold), mark, sew & cut as before.